



Pregnancy Test

Instructions: Midstream


- It is recommended to perform the test on the first time you urinate in the morning, as that is when hCG is most concentrated.
- Check that the test is completely sealed then open the pouch by tearing along the notch.
- Remove the cap from the device, then hold the device by the handle with the exposed absorbent tip pointing downwards & toward the body.
- Begin to urinate. Place the absorbent tip into the urine stream for at least 5 seconds, making sure that the entire tip is wet. Do not urinate past the blue mark.
- Remove the device from the urine stream. Replace the cap over the wick.
- Place the device on a flat surface with the result window facing up.
- As the test develops, the control or C line will appear quickly; this indicates the test is valid. If no C line appears, discard test and begin again with a new test.
- Test results will appear within 3 to 5 minutes. A test with two lines indicates a positive or pregnant result. A test with only the control or C line indicates a negative or not pregnant result. A test with no lines or only a T line result is invalid. Discard and begin again with a new test.
- Test results should not be read after 5 minutes.

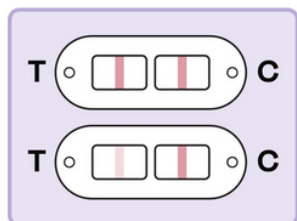

Positive Result

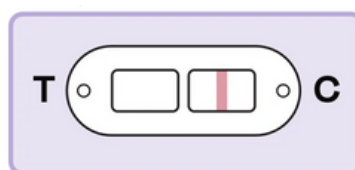


Negative Result

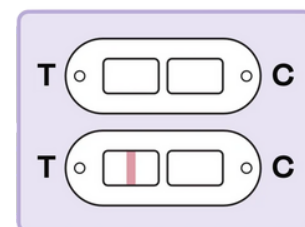


Invalid Result

Pregnancy Test Instructions: Cup


- It is recommended to perform the test on the first time you urinate in the morning, as that is when hCG is most concentrated.
- Check that the test is completely sealed then open the pouch by tearing along the notch.
- Urinate into a clean cup.
- Remove the cap from the device, then dip the into the urine. Do not allow the urine to reach the level of the blue mark.
- Replace the cap over the wick of test.
- Place the device on a flat surface with the result window facing up.
- As the test develops, the control or C line will appear quickly; this indicates the test is valid. If no C line appears, discard test and begin again with a new test.
- Test results will appear within 3 to 5 minutes. A test with two lines indicates a positive or pregnant result. A test with only the control or C line indicates a negative or not pregnant result. A test with no lines or only a T line result is invalid. Discard and begin again with a new test.
- Test results should not be read after 5 minutes.

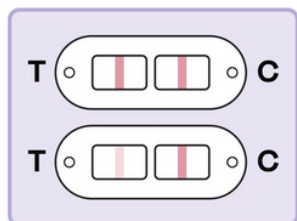

Positive Result

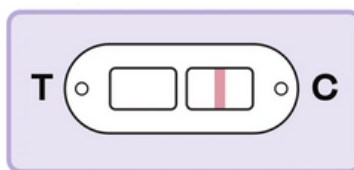

Negative Result

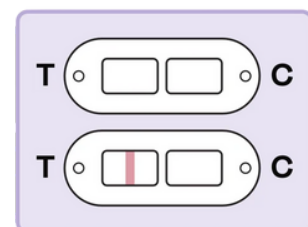

Invalid Result